

Kosher Circumcision Tray

Catalog # 00933


Convenient, Disposable, Sterile

Contains all components for required Jewish Kosher Circumcision procedure, according to "*Halacha*" – with Rabbinic approval.

Disposable instruments and tray eliminate risk of Cross- Infection and Nosocomial Infections.

The Tray Contains:

Waste bag with adhesive strip. The waste bag can be attached to the working surface for easy access during the procedure.

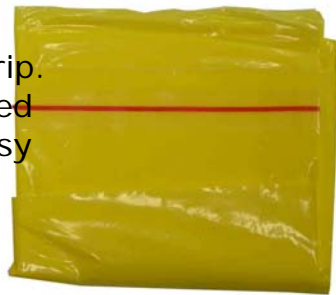

Three Benzalkonium Chloride **Towellettes** to **disinfect** the "*mohel*"'s hands and the infant's penis prior to the procedure.



Jumbo **Cotton Tipped Applicator** to administer Wine to the infant.



Europlas[®] Probe for easy separation of Prepuce from Glans Penis.

Stainless Steel Super Fine **Mosquito Clamp**



Clinically proven Stainless Steel extremely thin (0.7mm) **Mogen Circumcision Shield**. 1.2mm to 0.4mm convergent groove with blunt inner edges.

Mogen Circumcision Shield prevents side slipping, protecting Glans Penis.

Mogen Clamp ensures accurate visibility and gauging of Prepuce and Mucosa to be excised.





Europlas® Scalpel sharpened at both edges "*Cherev Pepiot*".

Sharp/Blunt Stainless Steel Sharp Tipped
Europlas® Scissors



Europlas® Forceps



Suction tube, transparent plastic material, made of two parts, one for the *Mohe!* (Circumciser) lips, the second for the penis. Between the two parts, an absolute filter that filters all types of microorganisms



Antiseptic Bandage, containing **silver Ions** (20%) that destroy any type of Microorganism.



The **silver** bandage is placed around the glans base, at the point of circumcision of the foreskin, the **absorbent compressed cotton** is rolled over it, the silicon paper is removed and the bandage is finally adhered by means of the **adhesive strip**.

4 **Gauze Sponges**.



Empty Bag to hold the Circumcised Foreskin

Pair of **latex gloves** (Optional)

Pollak Ltd. Hamelacha 20, Rosh Ha'ayin, Afek, 48091 Israel.
Tel: 972 3 9022885 Fax: 972 3 9024531 euroband@hotmail.com
www.euroband.com

EUROBAND

